
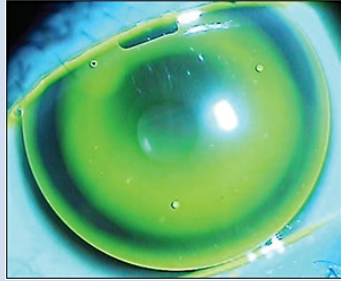
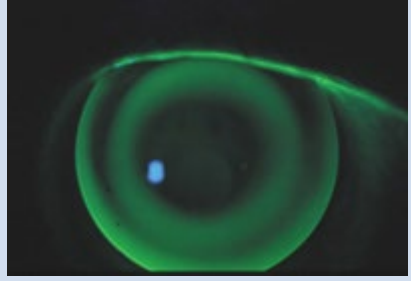


	Optimal	Steep	Flat
Fluorescein Pattern			
Centration	Good	Good may sit low	May sit laterally/inferiorly
Movement	Adequate and smooth	Minimal	Excessive
Central zone	Central bearing 3-5 mm Ø	Small area of bearing	Large area of bearing
Reverse zone	0.5-1 mm annulus of tear reservoir	Wide or barely visible with bubbles present	Optimal to wide annulus of tear reservoir
Alignment zone	1.0 – 1.5 mm annulus of mid peripheral bearing	Optimal to narrow	Narrow
Peripheral zone	Clearance approx. 0.4 mm wide	Narrow peripheral clearance	Wide peripheral clearance